

Detection of Nucleoli Using ImageJ

Nico Hochberger

November 2014

Contents

1	Introduction	5
1.1	Preface	5
1.2	Objective	5
1.3	Motivation	5
2	Design and Implementation	6
2.1	Requirements	6
3	User's Manual	9
3.1	System Requirements	9
3.2	Starting the Application	9
3.3	Configuration	9
3.3.1	General Parameters	9
3.4	General Picture Analysis Settings	10
3.4.1	Improved Image Detection Parameters	12
3.5	Structure of Files and Folders	12
4	Image Analysis	13
4.1	Target Detection Procedure	13
4.1.1	Determining Regions of Interest	13
4.1.2	Finding Nucleoli	13
4.2	Thresholding	13
4.3	Particle Analysis	13
5	Conclusion and Prospect	15
5.1	Conclusion	15
5.2	Prospect	15
	Bibliography	16

List of Figures

1	Example analysis performed with CellProfiler	5
2	Example of the image containing the results	8
3	Structure of directories and files created in the application's folder and the source folder	13
4	Schematic of the application's execution process	14
5	Example images of channel 1 and 3	14

List of Tables

1	Minimum system requirements of Java 1.7 on Windows 64bit	9
---	--	---

1 Introduction

1.1 Preface

1.2 Objective

1.3 Motivation

Currently nucleoli are detected using an application called CellProfiler¹. While this application yields reliable results, it also takes pretty long to complete the analysis. Runtimes up to 45 seconds are common.

Due to the fact that CellProfiler is a very general approach, applicable to a large variety of tasks related to detecting nuclei and nucleoli, the results of its analysis have to be checked manually to reduce the amount of false-positives.

- TODO -

Consequently, this leads to a more specialized way of analyzing the cells, which does not do all the analysis performed by CellProfiler, but on the other hand is much faster and thus helps to prevent cells dying before the analysis is finished.

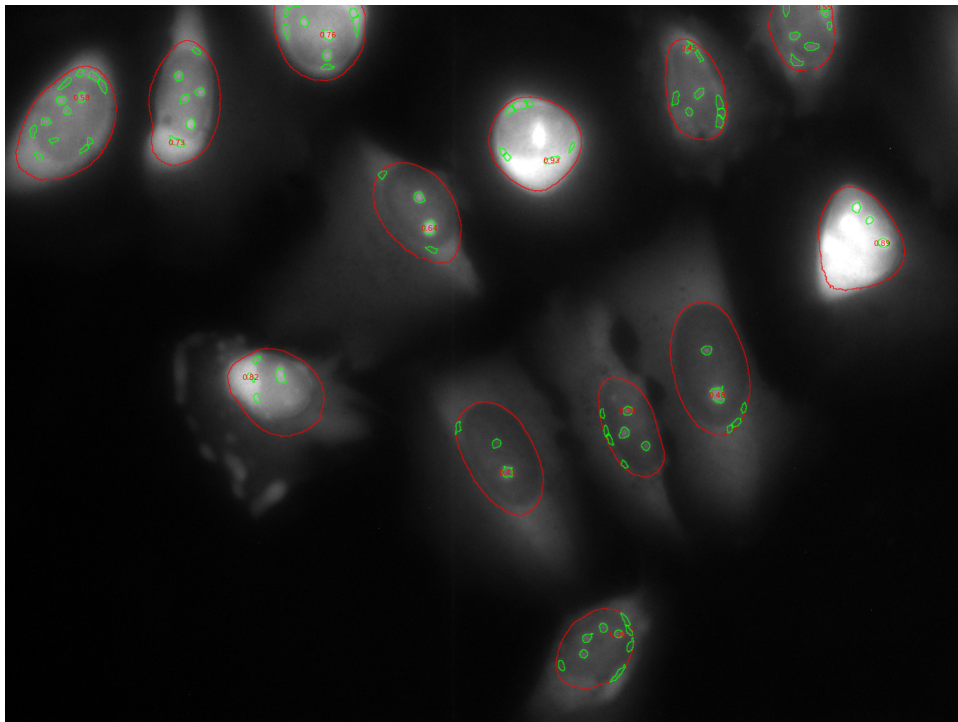


Figure 1: Example analysis performed with CellProfiler

¹<http://www.cellprofiler.org>

2 Design and Implementation

2.1 Requirements

The application has to meet the following requirements:

- **Reliable nucleoli detection:** Stable and reliable detection of nucleoli is the main purpose of the application. Hence, it is supposed to detect at least 75% of the nucleoli CellProfiler can detect. This includes a certain degree of stability concerning fuzzy pictures or pictures with inequally distributed or diffuse brightness.
- **Fast analysis:** As the application is tailored to this single task, it is expected to detect a suitable amount of nucleoli in only a small fraction of the time required by CellProfiler. The topmost time that the application may require to complete the analysis of one image is five seconds.
- **Fallback in case of empty nuclei:** Each nucleus is expected to contain at least one nucleolus. Yet, this expectation cannot always be fulfilled due to potentially damaged nuclei, fuzzy images, or other reasons. In this case, the center of the nucleus has to be provided as fallback target.
- **Visualizer:** In order to quickly check the results directly after running the analysis and to provide a way to quickly present the results to a potential audience, the application has to provide the possibility to be configured so that it shows the results as an image. This image has to contain all detected nucleoli targets, fallback targets and the regions of interest, e.g. the nuclei.
- **Versatility:** Since the appearance of different specimen can vary in various ways, all analysis parameters have to be configurable. Among others, this includes the minimum and maximum sizes of nuclei and nucleoli. The configuration is supposed to be achieved via an understandable, interchangeable, text-based file².
- **Statistics:** To determine the most suitable parameters for different kinds of specimen, another feature may be configured. This statistics feature has to include:
 - The amount of detected nuclei
 - The amount of detected nucleoli
 - Nuclei to nucleoli ratio as percentage

²Configurable parameters are explained in detail in the User's Manual section

- The distance of each detected nucleolus to the center point of the containing nuclei and the average distance in pixels
 - The area of each detected nucleus and the average area in square pixels
 - The area of each detected nucleolus and the average area in square pixels
- **Serialization of the results:** All results have to be stored in their accordant files in a subfolder *results* of the folder containing the original data. The name of each result file has to contain a the timestamp indicating the application's execution in the format `<year>_<month>_<day>_<hour>_<minute>_<second> .`
E.g. `targets_2014_11_20_14_50_35.txt` .

In the following, the accordant formats and files are described.

- **Targets:** Real targets and fallback targets sre to be saved in one txt-file named `targets_<timestamp>.txt` in the following format:

```

1 # nucleoli targets
2 <target number> : [<x-coord>, <y-coord>]
3 ...
4 # targets in center of empty nuclei
5 <target number> : [<x-coord>, <y-coord>]
6 ...

```

Listing 1: Format of results txt-file

Example:

```

1 # nucleoli targets
2 1 : [468, 43]
3 2 : [1183, 14]
4 # targets in center of empty nuclei
5 3 : [87, 174]
6 4 : [769, 198]

```

Listing 2: Example of results file

- **Statistics:** The statistics as mentioned above have to be stored in a txt-file named `statistics_<timestamp>.txt` . An example of the statistics file can be found in the appendix.
- **Result image:** The image as it would be displayed by the visualizer has to be stored to a file named `targets_<timestamp>.<original image filetype>` . See Figure 2 for an example the image.

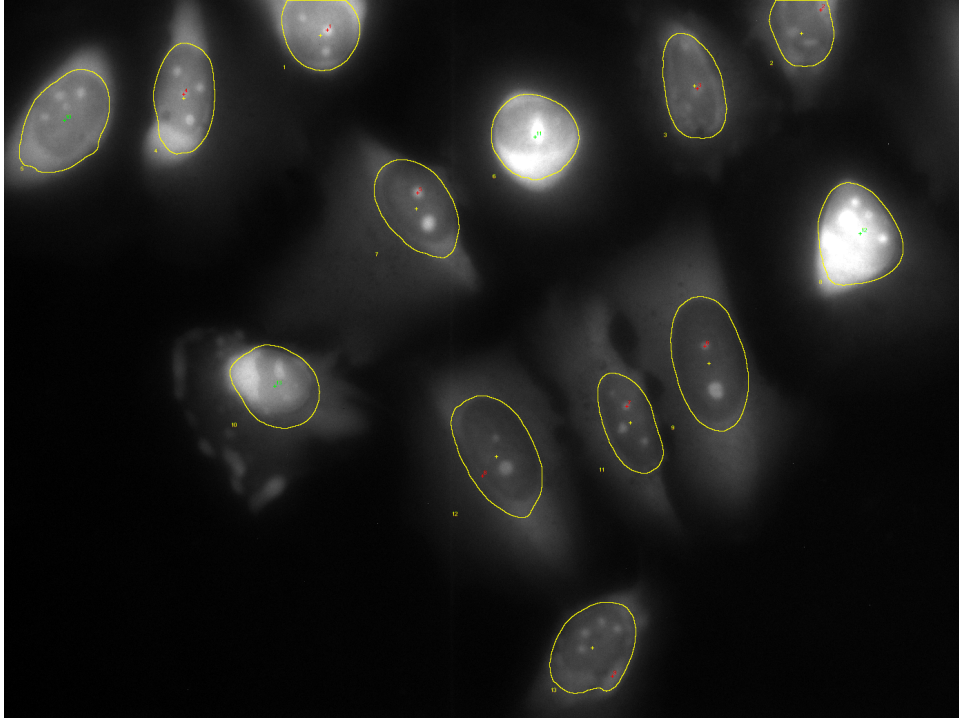


Figure 2: Example of the image containing the results

- **Quick and easy deployment:** The application is supposed to be deployed and run as easily as possible. Since the environment and the machines this will happen on may vary drastically, a versatile and independent way to do this is required. The Java framework delivers this possibility. Hence, the application will be developed using Java.

3 User's Manual

3.1 System Requirements

The application requires the Java Runtime Environment 1.7 or higher. For Microsoft Windows 64bit operating system this results in the following minimum system requirements ³:

Memory:	128MB
Processor:	Pentium 2 @ 266MHz
Disc space:	124 MB + 2 MB

Table 1: Minimum system requirements of Java 1.7 on Windows 64bit

As it is good practice and vital to the security of any system, the Java framework should be updated on a regular basis. Note that the requirements may change with any update.

3.2 Starting the Application

Starting the application is accomplished using the command line interface via Java's `java -jar` or `javaw -jar` command⁴. Furthermore, the application expects the configuration file to be passed as single parameter.

Considering a scenario in which the configuration file (see 3.3) is stored as `configuration.txt` in the same folder as the `NuFi.jar`, the following command will properly start the application:

```
java -jar NuFi.jar configuration.txt.
```

3.3 Configuration

Configuration of the application is performed using a simple text file containing the information listed in the following. Note that a missing parameter will result in a fatal error, preventing the application from properly working.

A complete example of the configuration file can be found in the appendix.

3.3.1 General Parameters

This section contains information on filetypes, source folders and naming convention of source files.

³Further details on the system requirements can be found at: <https://docs.oracle.com/javase/7/docs/webnotes/install/>

⁴See <https://docs.oracle.com/javase/7/docs/technotes/tools/windows/java.html> for further details.

- **Source folder:** The folder containing the source files (e.g. image files) on which the analysis will be based. The source folders are expected to be provided as absolut path or relative to `NuFi.jar`.

Note that the application expects exactly three files of the type defined as channel filetype. If there are more or less than three files of that type, a fatal error will occur.

Key:

`source.folder`

Examples:

`source.folder = C:/example/images`

`source.folder = ../../example/images`

`source.folder = example/images`

- **Used channels:** Each specimen is photographed three times using different colorization. This results in three different images referred to as different channels. Channels 1 and 3 are used in the process of detecting targets, thus, the correct specification and order of the files is vital for image analysis. Position one determines channel 1, position 3 determines channel 3. The source folder is scanned for files containing these channels determination. If these are not found, the application will terminate with a fatal error.

Key:

`used.channels`

Examples:

`used.channels = Kanal1, Kanal2, Kanal3`

`used.channels = channel1, channel2, channel3`

- **Channel filetype:** Though using the png-file format is advised, the application is built to be able to analyze various image formats. If it is necessary to analyze images that are not provided in the png-format, the accordant format can be provided using this key.

Key:

`channel.filetype`

Examples:

`channel.filetype = png`

`channel.filetype = jpg`

3.4 General Picture Analysis Settings

This section contains information of the expected size and shape of nuclei and nucleoli, as well as the width of the in depth-analysis.

- **Size settings:** Since the size of the nuclei and nucleoli in the specimen may vary and since the magnification of the microscope used in taking the pictures may also change, the size in the photographs may need to be adapted. This can be achieved by giving an average size in square-pixels and the minimum and maximum sizes as multiples of that size.

Keys:

```
nucleus.average  
nucleus.min.factor  
nucleus.max.factor  
nucleolus.average  
nucleolus.min.factor  
nucleolus.max.factor
```

Examples:

```
nucleus.average = 10000  
nucleus.min.factor = 0.5  
nucleus.max.factor = 1.5  
nucleolus.average = 100  
nucleolus.min.factor = 0.75  
nucleolus.max.factor = 1.8
```

- **Minimum circularity:** This parameter pays respect to the variation of the objects in their shape. A value of 1 indicates a perfect circle, while a 0 indicates that there are no requirements to the shape of the object. Since nuclei typically are elliptical in shape, a value of approximately 0.6 is advised. Nucleoli are typically close to a perfect circle. Thus, a value of 0.9 includes most of them while it excludes a variety of anomalies and prevents false-positives on nuclei borders.

Keys:

```
nucleus.min.circularity  
nucleolus.min.circularity
```

Examples:

```
nucleus.min.circularity = 0.6  
nucleolus.min.circularity = 0.9
```

- **In-depth width:** This parameter defines the variation in thresholding used during in-depth analysis.

Key:

```
indepth.range
```

Examples:

```
indepth.range = 25  
indepth.range = 10
```

3.4.1 Improved Image Detection Parameters

In the standard setting the application will use the improved image detection algorithm. In this section, the parameters required for this algorithm can be defined. These include setting for brightness correction (`<object>.background.blur` and `<object>.thresholding.blur`) and the value used for correcting discrepancies between channel 3 colorization and the acutal size of nuclei (`nucleus.boundary.width`).

- **Brightness correction:**

Keys:

```
nucleus.background.blur
nucleus.thresholding.blur
nucleus.boundary.width
nucleolus.background.blur
nucleolus.thresholding.blur
```

Examples:

```
nucleus.background.blur = 100
nucleus.thresholding.blur = 3
nucleus.boundary.width = 5
nucleolus.background.blur = 10
nucleolus.thresholding.blur = 1
```

3.5 Structure of Files and Folders

The results of the image analysis will be stored into a directory `results` which is created as a subfolder of the directory defined as source folder. This directory will contain the targets file, the result image and, in case statistics are enabled, the statistics file. See figure 3.

The application will neither delete nor overwrite files created by a previous execution since each file contains a timestamp as defined in 2.1.

Besides this, the application will create a directory `logs` on the same level as `NuFi.bat` is located. This directory contains the log-file of the application itself. If any unexpected behaviour is encountered while using the application, refer to this file since it might contain helpful information. See figure 3.

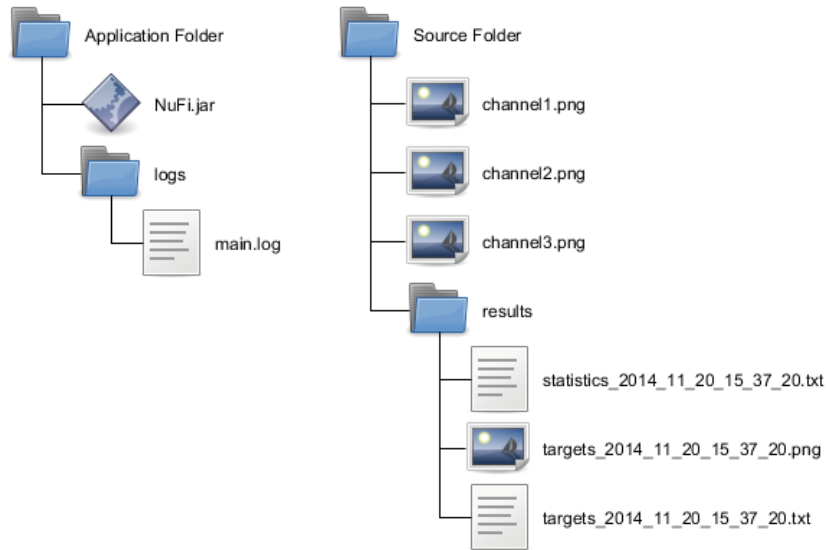


Figure 3: Structure of directories and files created in the application's folder and the source folder

4 Image Analysis

In this section the procedure of how the application processes and analyzes images is explained. Details such as how the application checks its configuration are omitted in favor of more emphasis on image analysis.

4.1 Target Detection Procedure

The target detection process can be basically divided into two parts. At first, the regions of interest are detected using channel 3. These regions are then transferred to channel 1 in order to only analyze those regions that may actually contain nucleoli.

In the following, these two parts are explained step-by-step using the example images in figure 5.

4.1.1 Determining Regions of Interest

4.1.2 Finding Nucleoli

4.2 Thresholding

4.3 Particle Analysis


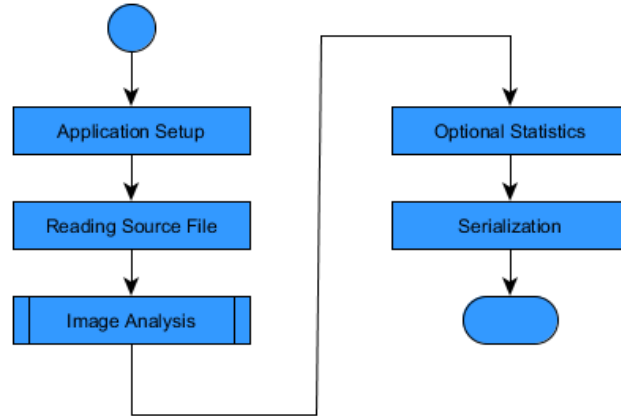


Figure 4: Schematic of the application's execution process

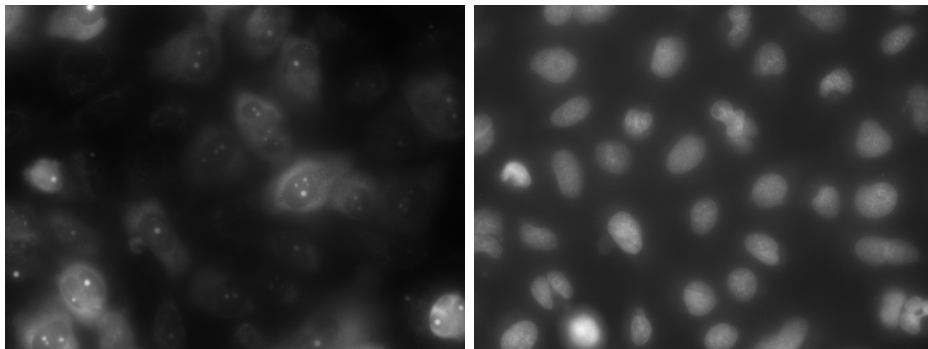

(a) Channel 1 example image

(b) Channel 3 example image

Figure 5: Example images of channel 1 and 3

5 Conclusion and Prospect

5.1 Conclusion

5.2 Prospect

References